



Demers eFX Ambulance

ALL-ELECTRIC TYPE III AMBULANCE

- Ambulance designed by and for First Responders
- Built on Lion5 all-electric chassis & designed to last 10 years
- Purpose-built ambulance ergonomically designed for EMS safety, comfort, & optimization
- Patient Compartment with state-of-the-art clinical design and 20%* larger to optimize crew comfort & mobility while providing optimal access to patients
- FX Action Seats give First Responders freedom of movement & accessibility to patients and equipment while belted

** Based on comparison to 164" ambulance module size.*

REVOLUTIONARY FX ACTION SEATS

FX ACTION SEATS

- Unique ergonomic design
- Maximizes mobility optimizing patient care while belted
- 10x more opportunities to remain safely buckled while caring for patient
- Stabilized seat to help perform critical tasks
- New design enables Medic to maneuver using legs, releasing hands to care for patients
- Installed near stretcher to prevent unbuckling
- Belt & occupant sensor for driver awareness

***patents pending*



Picture of ambulance may vary

eFXambulance.com

DEMERS IS DRIVEN TO IMPROVE CREW EXPERIENCE & LIFE SAVING OPERATIONS

Demers Ambulances USA Inc. • 1170 Production Drive • Van Wert, OH 45891 • Toll Free (800) 363-7591

A REIMAGINED AMBULANCE FEATURING A NEVER-BEFORE-SEEN CLINICAL DESIGN.



MEDIC FOCUSED PATIENT COMPARTMENT FEATURES

MODULE SIZE

- 184" long x 96" wide with 20% more medical space than the 164" model

PATIENT COMPARTMENT

- New and innovative, state-of-the-art design
- Focusing on clinical and ergonomic requirements
- Keeping safety & comfort as top priorities
- 72" interior headroom
- Less noise and vibration

WORK STATION WITH SWIVELING DRAWERS

- Promotes natural workflow, everything within reach of FX Action Seats
- Adjustable positioning, equipment remains secured with electrical release system
- Pivoting writing surface above drawer

ONE CLICK 4-POINT SEATBELT HARNESS

- Convenient to use with only one buckle to click
- 4-point harness holds you in securely
- Allows the freedom to stand while remaining buckled†

POWERLOAD EXTENDER

- Stryker PowerLOAD with Demers automated positioning system
- Primary care position & loading position near the rear doors
- Raises cot 3" to improve ergonomics & comfort while treating patients

CEILING CONSOLE

- O₂ port & IV bag hooks
- Multi-directional cooling vents above stretcher
- Convenient grab rails

LINEAR LIGHTING

- Reduces stress for patient with comforting and soothing lighting options
- (3) zone adjustable lighting (4 intensities) and (2) temperature color options

PASSTHRU WINDOW

- Ultrawide, electric window optimizes visual communication amongst crew members
- Optimizes privacy & reduces noise (seal)
- Independent climate control enhances comfort for cab and module occupants



†Standing in a moving vehicle is unsafe.



eFX EXTERIOR FEATURES

MMC MULTIPLEX ELECTRICAL SYSTEM

- System is proven to reduce downtime & lower repair costs
- Supports remote diagnostics & preventative maintenance
- Does not require certified electrical technical to perform maintenance

AERODYNAMIC CAB RISER

- Reduces energy consumption & increased visibility with emergency lighting
***patent pending*
- Integrated HVAC with easy servicing, no special tools required to access the system
- HVAC designed for efficiency & durability with optimal drainage

ADAPTIVE SUSPENSION

- Smooth, quiet, comfortable ride
- 3.2" kneeling capacity lowers step height to help load & unload patients
- Standard front & rear axle with independent controls

140" WHEELBASE

- Tight turning radius of 22.7 ft
- Easily navigate tight urban streets
- Optimal maneuverability with shortest ambulance wheelbase in class

FENDER

- Highly durable, UV resistant & easily replaceable with automotive fasteners

CHARGING

- Charging port curbside (front corner):
 - SAE J1772 Combo (CCS1) Connector
 - Type: DC Fast Charger (level 3)
 - Maximum charging power 100 kW
 - Recharge time up to 100 kW:
 - From 40% to 80% - 47 minutes
 - From 80% to 90% - 15.5 minutes
 - From 90% to 100% - 91 minutes (1.5 hours)
 - Output voltage: 800V DC
 - Lithium-Ion Battery Pack 800V





eFX CAB FEATURES

ERGONOMICALLY DESIGNED FOR EMS COMFORT, HEALTH & SAFETY

- Spacious cab with improved reclining space
- Optimal leg & shoulder room
- Category-leading front & side visibility
- Cup holders & storage space for personal belongings
- Easy access to electronic interfaces & charging ports

15" TOUCHSCREEN PANEL


- More functions, quicker access with the touch of a button
- Includes emergency functions, lights, sirens, climate control, 360 cameras & more
- Mirror cell phones, computers or tablets onto the display

STEERING WHEEL CONTROLS

- First ambulance steering wheel to incorporate features designed explicitly for EMS transport
- Emergency lights and switches on the steering wheel to help crew remain focused driving




WEBSITE

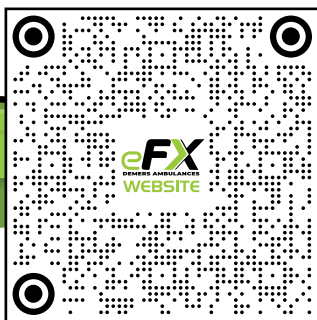

▲ DETAILS ▲

VIDEO

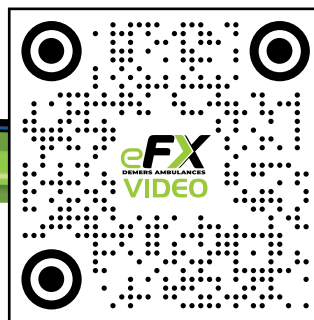

▲ WATCH NOW ▲

MORE RESOURCES

**Patent Pending. Publ. No. US 20240300402A1 (Aerodynamics Cab).
 Patent Pending. Publ. No. US 2023166641-A1 (Multi-Pivot Seat Base Assembly).
 Patent Pending. Publ. No. US 20240190303-A1 (Safety Seat for Emergency Vehicles).

eFXambulance.com

Demers Ambulances USA Inc. • 1170 Production Drive • Van Wert, OH 45891 • Toll Free (800) 363-7591